

## Agar Dilution Method for Detection of Inducible Clindamycin Resistance in *Staphylococcus* spp.<sup>∇</sup>

Clarence J. Fernandes,<sup>1</sup> Matthew V. N. O'Sullivan,<sup>1,2\*</sup> Yongwei Cai,<sup>2,3</sup> Fanrong Kong,<sup>2</sup>  
Xianyu Zeng,<sup>2,4</sup> Gwendolyn L. Gilbert,<sup>2</sup> and George Kotsiou<sup>1</sup>

Department of Microbiology, Royal North Shore Hospital, Sydney, New South Wales, Australia<sup>1</sup>; Centre for Infectious Diseases and Microbiology—Public Health, Institute of Clinical Pathology and Medical Research, Westmead, New South Wales, Australia<sup>2</sup>; Department of Dermatology, Hangzhou Third People's Hospital, Hangzhou, Zhejiang Province, People's Republic of China<sup>3</sup>; and Department of Dermatology, Wuhan First Hospital, Wuhan, Hubei Province, People's Republic of China<sup>4</sup>

Received 9 June 2007/Returned for modification 27 August 2007/Accepted 8 October 2007

**We describe the development and validation of an agar dilution method for the detection of inducible clindamycin resistance by using 227 previously characterized erythromycin-resistant, clindamycin-susceptible *Staphylococcus* sp. isolates. Mueller-Hinton agar with defibrinated horse blood containing a range of erythromycin concentrations (1 to 8 mg/liter) combined with clindamycin at 0.5 mg/liter was used to determine the optimal concentration that produced growth of inducible isolates while inhibiting that of isolates without the inducible phenotype. A concentration of clindamycin of 0.5 mg/liter with erythromycin at 1 mg/liter was the optimal combination for detection of inducible resistance and resulted in a sensitivity of 100% (95% confidence interval [CI], 97.9 to 100) and a specificity of 100% (95% CI, 93.0 to 100). Attention must be paid to ensuring that a sufficient inoculum has been used, since an inoculum below the standard 10<sup>7</sup> bacteria/ml may result in false-negative results. This method has been incorporated into routine use in our laboratory.**

Macrolide resistance in *Staphylococcus* spp. is generally attributable to one of two mechanisms (7). One is active efflux due to a pump encoded by the *msrA* gene, which confers resistance to macrolides and streptogramin antibiotics but not lincosamides such as clindamycin (MS phenotype). The second is ribosomal methylation mediated by enzymes encoded by one of a variety of *erm* genes. This mechanism results in resistance to macrolides, lincosamides, and streptogramin B antibiotics (MLS<sub>B</sub> phenotype). The MLS<sub>B</sub> phenotype is inducible (iMLS<sub>B</sub>) or constitutive (cMLS<sub>B</sub>). Induction of resistance occurs in the presence of macrolide antibiotics but not lincosamides. Thus, isolates exhibiting iMLS<sub>B</sub> resistance appear susceptible to clindamycin by routine methods such as broth microdilution, disk diffusion, and agar dilution. However, treatment failures have occurred with such isolates when clindamycin is used, because of selection of mutations during therapy, leading to cMLS<sub>B</sub> resistance (4, 8, 12).

It is therefore important for microbiology laboratories to show that an isolate that is erythromycin resistant but apparently clindamycin susceptible on routine testing is of the MS phenotype before reporting it as clindamycin susceptible. This is most commonly done by using the disk diffusion test, in which a 2- $\mu$ g clindamycin disk is placed next to a 15- $\mu$ g erythromycin disk (D-zone test). Flattening of the clindamycin zone adjacent to the erythromycin disk (positive D-zone test) indicates the presence of iMLS<sub>B</sub> resistance, while a circular zone (negative D-zone test) indicates MS resistance. This method has been shown to have high sensitivity and specificity com-

pared with genotypic analysis and has been incorporated into the Clinical and Laboratory Standards Institute (CLSI) method (3), but false-negative tests may occur if the disk separation distance is too wide (5, 9).

For laboratories that use agar dilution for routine antimicrobial susceptibility testing, disk diffusion methods for detection of inducible clindamycin resistance are inconvenient. An agar dilution method has the potential benefit of allowing large numbers of isolates to be simultaneously tested. We describe the development and validation of an agar dilution method for detection of inducible clindamycin resistance that has been incorporated into routine use at our laboratory.

### MATERIALS AND METHODS

**Collection and characterization of isolates.** Consecutive, nonduplicate *Staphylococcus* sp. isolates which were clindamycin susceptible (MIC,  $\leq$ 0.5 mg/liter) but erythromycin intermediate (MIC, 1 to 4 mg/liter) or resistant (MIC,  $\geq$ 8 mg/liter) by broth microdilution were collected, identified, and stored as outlined in a previous study (9). The D-zone test was performed according to the CLSI method with Mueller-Hinton agar as described elsewhere (3, 9), and genes conferring macrolide resistance, namely, *ermA*, *ermB*, *ermC*, *ermTR*, and *msrA*, were identified with a multiplex PCR-based reverse line blot assay as previously described (9). Three isolates used in the previous study were nonviable, and a fourth, with iMLS<sub>B</sub> resistance by D-test and no erythromycin resistance gene detected, was excluded from the analysis. Of the 227 remaining isolates, 176 had the iMLS<sub>B</sub> phenotype (100 methicillin-resistant *Staphylococcus aureus*, 56 methicillin-sensitive *S. aureus* [MSSA], and 20 coagulase-negative *Staphylococcus* sp. strains) while 51 had the MS phenotype (3 methicillin-resistant *S. aureus*, 3 MSSA, and 45 coagulase-negative *Staphylococcus* sp. strains). All isolates with the MS phenotype harbored *msrA* alone, while all those with the iMLS<sub>B</sub> phenotype harbored either *ermA* (92 isolates) or *ermC* (84 isolates). Two with *ermA* also harbored *msrA* and had the iMLS<sub>B</sub> phenotype. All isolates were erythromycin resistant by broth microdilution except one MSSA strain which was D-zone test positive and carried the *ermA* gene.

**Agar dilution.** Mueller-Hinton agar (BBL, Sparks, MD) with 3.3% defibrinated horse blood (MHA-HB), containing 1, 2, 4, or 8 mg/liter erythromycin and 0.5 mg/liter clindamycin (Upjohn Laboratories, Kalamazoo, MI) was prepared. In addition, agar plates with 0.5 mg/liter clindamycin alone or with 1

\* Corresponding author. Mailing address: Centre for Infectious Diseases and Microbiology, Level 3, ICPMR, Westmead Hospital, P.O. Box 533, Wentworthville, 2145 New South Wales, Australia. Phone: (612) 98455973. Fax: (612) 98938659. E-mail: matthew.osullivan@swahs.health.nsw.gov.au.

<sup>∇</sup> Published ahead of print on 17 October 2007.

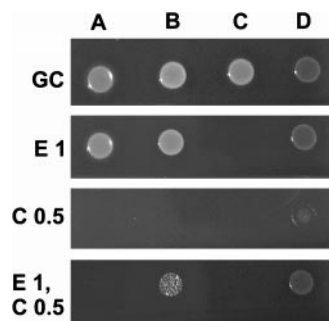


FIG. 1. Agar dilution method for detection of inducible clindamycin resistance. Columns: A, MS strain (ATCC BAA-976); B, iMLS<sub>B</sub> strain (ATCC BAA-976); C, macrolide-sensitive strain; D, cMLS<sub>B</sub> strain. Rows: GC, growth control (no antibiotics); E 1, erythromycin at 1 mg/liter; C 0.5, clindamycin at 0.5 mg/liter; E 1, C 0.5, erythromycin at 1 mg/liter and clindamycin at 0.5 mg/liter.

mg/liter erythromycin alone and agar plates without antibiotics were prepared, the latter two serving as growth controls. Mueller-Hinton agar supplemented with blood is used in our laboratory for all agar dilution susceptibility testing so that all organisms can be tested on the same plates. Additionally, we find that the opaque red background facilitates the detection of scanty growth. MHA-HB was poured into 100-mm-square plates to a depth of 3 mm. Once prepared, plates were stored at 4°C and used within 7 days. Approximately four colonies of an 18- to 20-h-old subculture were inoculated into 4 ml of Trypticase soy broth (BBL, Sparks, MD). After 4 h of incubation at 37°C, this suspension was adjusted to a concentration equivalent to a 0.5 McFarland standard. A 100-μl volume of this suspension was then added to 900 μl of broth contained in a 64-well seed block to produce a concentration of approximately 1 × 10<sup>7</sup> bacteria/ml. Plates were then inoculated with an automated multipoint inoculator (Denley-Tech, Billingshurst, United Kingdom) which delivers a volume of approximately 1 μl per spot, resulting in an estimated final inoculum of 1 × 10<sup>4</sup> bacteria. Plates were incubated for 18 h at 35°C under atmospheric conditions. Growth was recorded if at least one colony was evident at the inoculation site after careful visual inspection.

**Test interpretation.** A test was deemed to be positive (i.e., inducible resistance was present) if there was any visible growth on the erythromycin-only and combined plates but not on the clindamycin-only plate. A test was deemed to be negative (i.e., MS resistance phenotype) if growth was found on the erythromycin-only plate but not on the combined or clindamycin-only plate (Fig. 1). Genotyping was used as the reference for comparison, where a negative result was indicated by an *msrA* genotype and no *erm* gene detected, while a positive result was an *ermA* and/or *ermC* genotype regardless of whether *msrA* was present.

**Quality control.** Negative and positive control strains—ATCC BAA-976 (MS phenotype, *msrA* gene positive) and BAA-977 (iMLS<sub>B</sub> phenotype, *ermA* gene positive)—and two in-house isolates—one erythromycin and clindamycin suscep-

tible and one erythromycin and clindamycin resistant (cMLS<sub>B</sub> phenotype)—were inoculated onto each plate.

**Inoculum study.** To determine the importance of an adequate inoculum, 11 *S. aureus* isolates and 1 *S. epidermidis* isolate with inducible clindamycin resistance (comprising 9 *ermA*-positive and 3 *ermC*-positive isolates), as well as the four control strains, were tested with dilutions of 10<sup>8</sup>, 10<sup>6</sup>, and 10<sup>5</sup> bacteria/ml in addition to the standard inoculum concentration of 10<sup>7</sup> bacteria/ml. These were inoculated onto the same set of MHA-HB plates, incubated, and inspected as described above. In addition, 10-μl loops were used to inoculate blood agar plates for colony counts to ensure the correct inoculum density.

**Statistics.** Confidence intervals were calculated by Confidence Interval Analysis for Windows (available at <http://www.medschool.soton.ac.uk/cia>).

RESULTS

Determination of the optimal erythromycin concentration.

The effects of increasing erythromycin concentrations on the sensitivity and specificity of inducible clindamycin resistance detection are shown in Table 1. The number of false-negative results increased at higher concentrations. All isolates with the MS phenotype were inhibited at all of the erythromycin concentrations tested. A combination of erythromycin at 1 mg/liter and clindamycin at 0.5 mg/liter was determined to be optimal for testing since this was the only combination that gave 100% sensitivity and specificity.

**Assay validation.** At antibiotic concentrations of 1 mg/liter erythromycin and 0.5 mg/liter clindamycin, all isolates with inducible resistance showed visible growth and all with the MS phenotype failed to grow, resulting in a sensitivity and specificity of 100% (Table 1). However, plates had to be inspected carefully for growth, which was often scanty, with as few as one colony visible on the inoculum for some isolates.

**Determination of inoculum effect.** There was a clear relationship between the inoculum density and the amount of growth visible on the plates. With erythromycin at 1 mg/liter and clindamycin at 0.5 mg/liter and a standard inoculum concentration of 1 × 10<sup>7</sup> bacteria/ml, all 12 isolates with inducible resistance showed visible growth, but in some instances this was as little as one colony. At 1 × 10<sup>8</sup> bacteria/ml, there was good growth of all of the isolates tested. Lower inoculum densities resulted in false-negative results—1 of 12 at an inoculum concentration of 1 × 10<sup>6</sup> bacteria/ml and 4 of 12 at 1 × 10<sup>5</sup> bacteria/ml.

TABLE 1. Effects of various concentrations of erythromycin on detection of inducible clindamycin resistance by the agar dilution method

Erythromycin concn (mg/liter) <sup>a</sup>	No. of strains with genotype of:		Agar dilution result	% Sensitivity <sup>b</sup>	% Specificity <sup>b</sup>
	<i>msrA</i> only (n = 51)	<i>ermA</i> and/or <i>ermC</i> (n = 176)			
1	51	0	Negative	100 (97.9–100)	100 (93.0–100)
	0	176	Positive		
2	51	3	Negative	98.3 (95.1–99.7)	100 (93.0–100)
	0	173	Positive		
4	51	5	Negative	97.2 (90.5–97.6)	100 (93.0–100)
	0	171	Positive		
8	51	9	Negative	94.9 (91.6–98.1)	100 (93.0–100)
	0	167	Positive		

<sup>a</sup> The clindamycin concentration was 0.5 mg/liter with an inoculum concentration of 10<sup>7</sup> bacteria/ml.

<sup>b</sup> 95% confidence intervals are in parentheses.

## DISCUSSION

Isolates of *Staphylococcus* spp. that are erythromycin resistant but clindamycin susceptible should not be reported as clindamycin susceptible unless iMLS<sub>B</sub> resistance has been excluded (4, 7, 11, 12). Currently, the only method routinely recommended for testing of iMLS<sub>B</sub> resistance is disk approximation testing (1–3). While this accurately distinguishes iMLS<sub>B</sub> resistance (5, 9), it is relatively labor intensive for laboratories which routinely use agar dilution methods for susceptibility testing or when large numbers of isolates are to be tested.

In this study, we have found that agar dilution testing with erythromycin at 1 mg/liter and clindamycin at 0.5 mg/liter is a sensitive and specific method for detection of iMLS<sub>B</sub> resistance in *Staphylococcus* spp. There was an obvious inoculum effect associated with iMLS<sub>B</sub> resistance, which is similar to the effect seen for some other antibiotic resistance mechanisms such as β-lactamase production (10, 13). As for all of the agar dilutions tests, the importance of the correct inoculum should not be underestimated (6). While using the standard inoculum concentration of  $1 \times 10^7$  bacteria/ml resulted in accurate test performance, growth was frequently scant (as little as one colony). Thus, it is important to ensure that a sufficiently dense inoculum is used. A higher inoculum concentration ( $1 \times 10^8$  bacteria/ml) could be used to produce more obvious growth, but this would not integrate easily into the standard CLSI susceptibility testing methods.

## ACKNOWLEDGMENT

Many thanks to Judy Kelly for laboratory assistance.

## REFERENCES

- Bell, S. M., B. J. Gatus, J. N. Pham, and D. L. Rafferty. 2006. Antibiotic susceptibility testing by the CDS method. [Online.] <http://web.med.unsw.edu.au/cdstest>. South Eastern Area Laboratory Services, Sydney, New South Wales, Australia.
- British Society For Antimicrobial Chemotherapy. April 2004. Testing for dissociated resistance in staphylococci. [Online.] [http://www.bsac.org.uk/\\_db/\\_documents/Testing\\_for\\_dissocit.pdf](http://www.bsac.org.uk/_db/_documents/Testing_for_dissocit.pdf). Accessed 24 October 2006.
- Clinical and Laboratory Standards Institute. 2006. Performance standards for antimicrobial disk susceptibility tests; approved standard—ninth edition. CLSI document M2-A9. Clinical and Laboratory Standards Institute, Wayne, PA.
- Drinkovic, D., E. R. Fuller, K. P. Shore, D. J. Holland, and R. Ellis-Pegler. 2001. Clindamycin treatment of *Staphylococcus aureus* expressing inducible clindamycin resistance. *J. Antimicrob. Chemother.* **48**:315–316.
- Fiebelkorn, K. R., S. A. Crawford, M. L. McElmeel, and J. H. Jorgensen. 2003. Practical disk diffusion method for detection of inducible clindamycin resistance in *Staphylococcus aureus* and coagulase-negative staphylococci. *J. Clin. Microbiol.* **41**:4740–4744.
- Jorgensen, J. H., and J. D. Turnidge. 2003. Susceptibility test methods: dilution and disk diffusion methods, p. 1111. In P. R. Murray, E. J. Baron, J. H. Jorgensen, M. A. Tenover, and R. H. Tenover (ed.), *Manual of clinical microbiology*, 8th ed., vol. 1. ASM Press, Washington DC.
- Leclercq, R. 2002. Mechanisms of resistance to macrolides and lincosamides: nature of the resistance elements and their clinical implications. *Clin. Infect. Dis.* **34**:482–492.
- Levin, T. P., B. Suh, P. Axelrod, A. L. Truant, and T. Fekete. 2005. Potential clindamycin resistance in clindamycin-susceptible, erythromycin-resistant *Staphylococcus aureus*: report of a clinical failure. *Antimicrob. Agents Chemother.* **49**:1222–1224.
- O'Sullivan, M. V. N., Y. Cai, F. Kong, X. Zeng, and G. L. Gilbert. 2006. The influence of disk separation distance on the accuracy of the disk approximation testing for inducible clindamycin resistance in *Staphylococcus* spp. *J. Clin. Microbiol.* **44**:4072–4076.
- Queenan, A. M., B. Foleno, C. Gownley, E. Wira, and K. Bush. 2004. Effects of inoculum and β-lactamase activity in AmpC- and extended-spectrum β-lactamase (ESBL)-producing *Escherichia coli* and *Klebsiella pneumoniae* clinical isolates tested by using NCCLS ESBL methodology. *J. Clin. Microbiol.* **42**:269–275.
- Rayner, C., and W. J. Munckhof. 2005. Antibiotics currently used in the treatment of infections caused by *Staphylococcus aureus*. *Intern. Med. J.* **35**(Suppl. 2):S3–S16.
- Siberry, G. K., T. Tekle, K. Carroll, and J. Dick. 2003. Failure of clindamycin treatment of methicillin-resistant *Staphylococcus aureus* expressing inducible clindamycin resistance in vitro. *Clin. Infect. Dis.* **37**:1257–1260.
- Syriopoulou, V. P., D. W. Scheifele, C. M. Sack, and A. L. Smith. 1979. Effect of inoculum size on the susceptibility of *Haemophilus influenzae* b to β-lactam antibiotics. *Antimicrob. Agents Chemother.* **16**:510–513.